

PerfectSeal[®]

Adhesive Wound Dressing


Healthy Solutions...

Hydrocolloid Wound Dressing


Product Features

Hydrocolloid is a material containing gel-form agents, such as carboxymethylcellulose(CMC) and gelatin. It is applied to wound management, support wound bed a moist environment, accelerate wound healing process.

1. CMC in the hydrocolloid convert in gel when in contact with wound exudates, accelerate wound healing.
2. Semi-transparent, better for observe wound healing process.
3. Elastic, comfortable, free movement, suitable to any body parts.
4. PU film surface, water-proof, reduce risk of roll up and bunching.

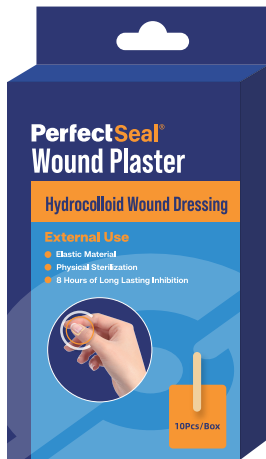


For the Hand

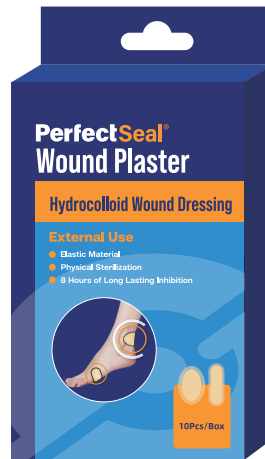


For the Heel



For the Toe

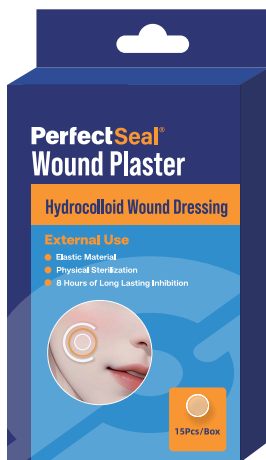

For the Face

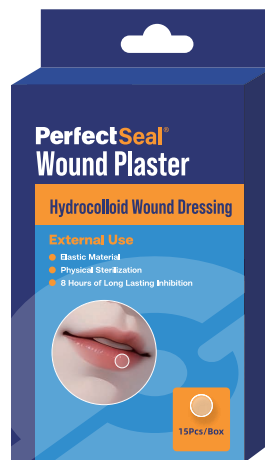

For the Lip



Product Parameter


Product	Item No.	Size(mm)	Package
	6501	55×37mm	10pcs/box, 100/case
	6502	30×30mm	10pcs/box, 100/case
	6503	20×65mm	10pcs/box, 100/case
	6504	20×60mm	10pcs/box, 100/case
	6505	42×62mm	10pcs/box, 100/case
	6506	30×30mm	10pcs/box, 100/case
	6507	44×69mm	10pcs/box, 100/case
	6508	80×120mm	10pcs/box, 100/case
	6509	150×180mm	10pcs/box, 100/case

■ Silicone Tape

Product Features

- 1.Gentle on the skin and painless to remove.
- 2.Reusable,does not stick to hair and leaves no residue, no remains.
- 3.Easy to use and hypoallergenic,with good adhesion.
- 4.Waterproof,but permeable.



Product	Item No.	Volume(mm)	Package
	6401	25mm×5m	12pcs/box, 20boxes/ctn
	6402	50mm×5m	6pcs/box, 20boxes/ctn
	6403	75mm×5m	4pcs/box, 20boxes/ctn

Other sizes are available based on customer's requests.

■ Silicone Foam Dressing

Product Features

Silione Foam is a foam dressing with silicone gel coating, designed to maintain a good healing environment for moderate to high exudate wounds such as pressure ulcers, leg and foot ulcers, traumatic wounds and other secondary healing wounds.

1. Gentle on the skin and painless to remove, with no residue.
2. "Vertically" absorbs wound mucosa exudate, reducing the risk of maceration around the wound.
3. Absorption capacity is great, this dressing keeps the wound under appropriate temperature and humidity according to basic "moist healing" principals.
4. Can be reused after inspection of the wound and will not adhere to the wound bed.
5. Does not damage the new epithelial cells when removed.
6. Helps to reduce scarring.
7. Waterproof, yet permeable. Showering is allowed while using this dressing.
8. Soft and comfortable, with good compression, works as a "cushion".

Product Parameter

Product	Type	Volume	Core Size
 	Silicone Non-Border	5 x 5cm	5 x 5cm
		7.5 x 7.5cm	7.5 x 7.5cm
		10 x 10cm	10 x 10cm
		15 x 15cm	15 x 15cm
	Silicone With-Border	5 x 5cm	2.5 x 2.5cm
		7.5 x 5cm	4.5 x 4.5cm
		8.5 x 5cm	5 x 2.5cm
		7.5 x 7.5cm	4.5 x 4.5cm
	Non-Border	5 x 5cm	-
		7.5 x 7.5cm	-
		10 x 10cm	-
	Self-Adhesive	5 x 5cm	2 x 2cm
		10 x 6cm	6 x 3cm
10 x 10cm		5 x 5cm	

Wound Dressing

Product Features

- 1.Provides a waterproof, sterile barrier to external contaminants including liquids, bacteria and viruses.
- 2.The dressing provides continuous visibility to the I.V. insertion site.
- 3.Adhesive is friendly on skin.
- 4.Dressing flexes with skin for greater patient comfort.
- 5.High breathable, the skin can breath freely.

Types



SecureFix plus



SecureFix



SecureFix+catheter securement device



SecureFix pad

Product Parameter

Name	Product	Size(cm)	Application Advice	Notes
SecureFix plus		8.5×11.5cm	CVC/PICC External Jugular Vein	Border design improve the adhesion,reduce edge-lift Dot coated border improve the Breathability Reduce the incidence of phlebitis caused by i.v.catheters Notes: without paper frame
		6.5×7cm	Arterial catheter,PICC,CVC,epidural catheter,PIV (peripheral intravenous catheter)	
		7×9cm	Arterial catheter,PICC,CVC,epidural catheter,PIV (peripheral intravenous catheter)	
		8.5×10.5cm	Peripheral and midline venous catheters,Central venous catheters (including subclavian jugular,femoral and PICCs),Implanted infusion devices(ports),Dialysis catheters,Pulmonary artery catheters,Arterial catheters	Non woven reinforced design improve the adhesion Dot coated border improve the Breathability Reduce the incidence of phlebitis caused by i.v.catheters Notes: without paper frame
		7×8.5cm	Peripheral lines	
		5×5.7cm	Peripheral lines Paediatric	
		7×9.5cm	For use on peripheral and central venous catheter fixation	
SecureFix+catheter securement device		10×12cm 9×3.5cm 6×7cm 9×3.5cm	Both secure the i.v.injection site and catheter	The fixation of the catheter reduces the medical risks,such as the retention needle falling off due to accidental pull-out.
SecureFix		6×7cm 7×9cm	Suitable for wing type or with pot retention needle fixation.	U notch designed basic type S peeling design on paper frame makes it easy to peel off the backing paper
		4.4×6cm	For Pediatric	Basic paper frame type S-peeling design on paper frame makes it easy to peel off the backing paper
		6×7cm	PIV(peripheral intravenous catheter)	
	10×12cm	CVC/CVP cervical veins,epidural catheters,central veins.		


Liquid Bandages

Product Features

Liquid bandage is a topical skin treatment for minor wounds that binds to the skin to form a protective polymeric layer that keeps dirt and germs out and moisture in. Both seal and protect minor cuts, scrapes, and burns. Liquid bandage typically forms a protective, waterproof layer over the wound to help it heal. At the same time, skin glue, also known as medical adhesive, is a type of surgical glue that is used to close small incisions or lacerations without the need for stitches.

1. Covers and Protects Wounds.
2. Kills 99.9% of germs to help prevent infection.
3. Protects cuts and scrapes.
4. Helps prevent and protect blisters.
5. Helps prevent the formation of calluses.



Product	Item No.	Size	Package
	9101	5g	1pcs/box, 24boxes/ctn



Product	Item No.	Size	Package
	9102	1g	2pcs/card, 12cards/box, 12boxes/ctn

PerfectSeal® DENSUN™

ZHEJIANG PERFECTSEAL NEW MATERIAL TECHNOLOGY CO.,LTD
NINGBO DENSUN MEDICAL TECHNOLOGY CO.,LTD

A2-1 of # 5 Plant , No.96 Mujin Road, Ningbo High-Tech Zone, Ningbo City, Zhejiang Province, China.

Tel: +86-574-8756 4800 / +86-574-87564717 E-mail: sales@perfectseal.com Http: www.perfectseal.com